

ORIPLAST[®]



Pleural Puncture

MEET PRECISE GOALS

Pleural Puncture

RECONSTRUCTION OF THE RESPIRATORY SYSTEM

Special Puncture Method

Utensils for pleural puncture are particularly vital in emergency medicine: They are used for immediate removal of e.g. fluid from the chest (pleural effusion) impairing the patient's breathing. In addition to such a therapeutic thoracocentesis, diagnostic pleural puncture is also important, which allows for the sterile extraction and analysis of the effusion.

Possible causes of pleural effusion are e.g. inflammatory changes, tumour diseases, traumatic injuries and heart diseases. The puncture is performed under local anaesthesia with the patient sitting or lying down.

The transparent, three-part syringe is designed for quick and easy connection to the cannula with a Luer-Lock connector. The syringe has a volume of 50 ml.

Syringe


Our pleural puncture set includes a 14G cannula. Thanks to the Luer attachment, this can be connected quickly and easily and is ready for use immediately.

Cannula

The adapter is suitable for the secure connection of a cannula.

Luer-Lock adapter



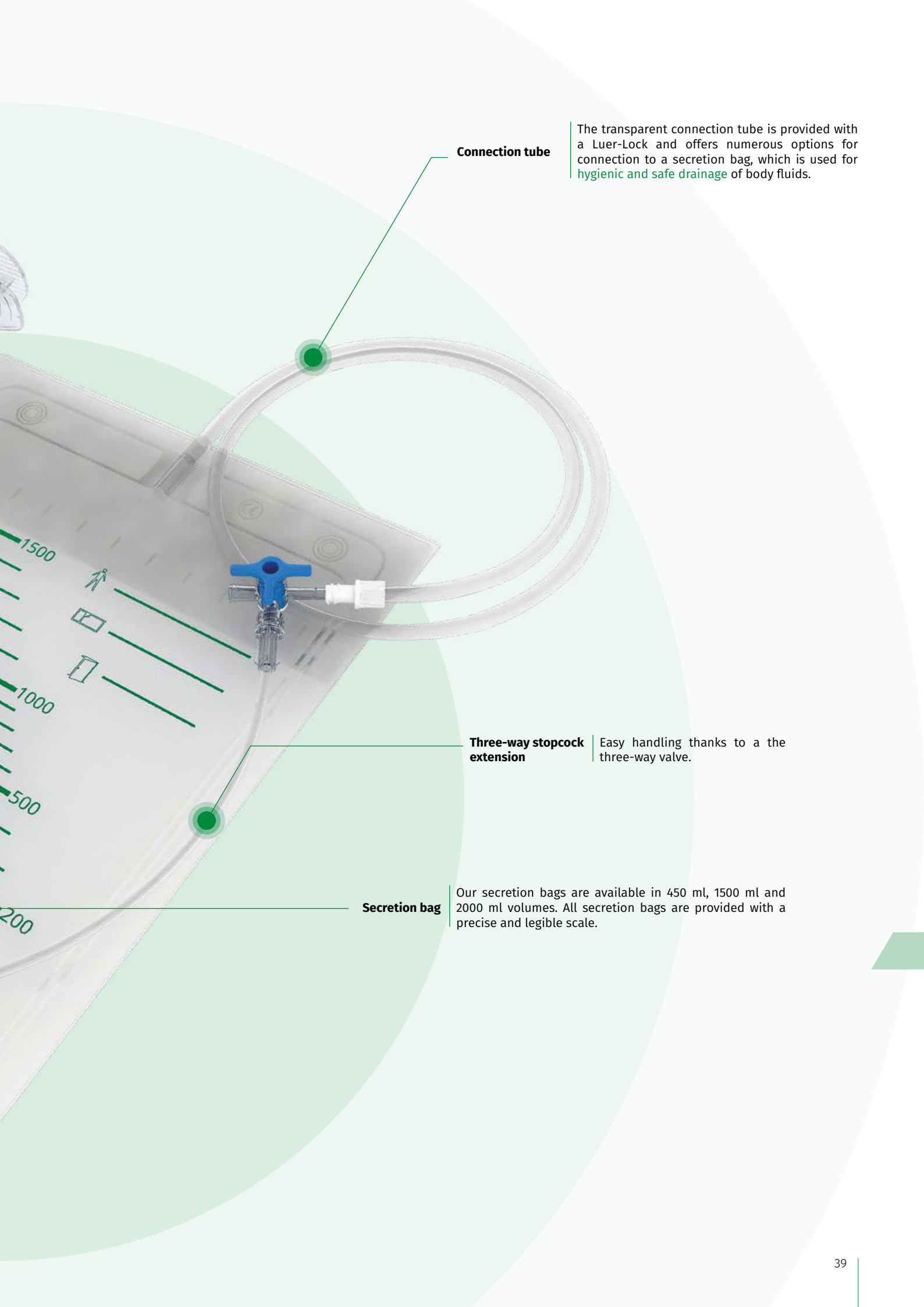

Connection tube

The transparent connection tube is provided with a Luer-Lock and offers numerous options for connection to a secretion bag, which is used for **hygienic and safe drainage** of body fluids.

Three-way stopcock extension

Easy handling thanks to a the three-way valve.

Secretion bag

Our secretion bags are available in 450 ml, 1500 ml and 2000 ml volumes. All secretion bags are provided with a precise and legible scale.

Pleural Puncture Set

This compact set is used in the diagnostic and therapeutic puncture of the intrapleural space e.g. in the presence of a pleural effusion.



Pleural Puncture Set

Secretion bag: tube length 900 mm, with non-return valve, LL pos, three-way stopcock (blue)

Syringe: 50 ml



Art.-No.	PCN	eClass	Volume (ml)	Bottom drain valve	Union cap	Cannula (G)	PU
VK00134	17484284	34-29-02-12	1500	-	-	14	50
VK00140	17484315	34-29-02-12	1500	-	•	14	50
VK00135	17484290	34-29-02-12	1500	-	-	-	50
VK00136	17484309	34-29-02-12	2000	•	-	14	50

Accessory

The exchange bag is suitable for collecting secretion fluids from the pleural cavity.



Secretion Bag

- tube length 900 mm



Art.-No.	PCN	eClass	Volume (ml)	Non-return valve	Bottom drain valve	Adapter	Protection cap	PU
VK00918	17488543	34-29-03-92	1500	•	-	LL neg	-	100
VK01070	17489502	34-29-03-05	1500	-	-	LL neg	-	100